



Mouse Anti-human Thr181 phosphorylated Tau Antibody Datasheet

Product Name: mAb anti-human Thr181 phosphorylated Tau	Clone No.: 9B4	Lot No.:
Catalogue No.: MO-M40099C	Quantity:	Expiration:

Description: Mouse monoclonal antibody against Thr181 phosphorylated Tau protein

Purification: Protein G affinity purified

Target Protein: Human Thr181 phosphorylated Tau

Immunogen: KLH conjugated with a short peptide [Ala-Pro-Lys-Thr(p)-Pro-Pro-Ser-Ser-Gly-Glu-Cys], corresponding to the amino acid sequence 177-187 on the phosphorylated Tau protein.

Fusion Myeloma: Sp2/0-Ag14

Specificity: This antibody reacts with the unconjugated Thr181 phosphorylated Tau peptide [Ala-Pro-Lys-Thr(p)-Pro-Pro-Ser-Ser-Gly-Glu-Cys]. Cross-reactivity with unphosphorylated peptide was undetectable or very low.

Species Reactivity: Human, mouse

Host / Isotype: Mouse/ IgG1, κ

Formulation: Lyophilized from a solution in 0.01M PBS pH7.2

Reconstitution: Double distilled water is recommended to adjust the final concentration to 1.00 mg/mL. Avoid repeated freeze and thaw cycle.


Storage: Store at -20°C

Research Area: Alzheimer’s disease and other neurodegenerative diseases

Background: Tau is a highly soluble microtubule associated protein found in high concentration in neuron of central nervous system. The main function of the protein is stabilizing the microtubule in nerve fibre (axon). There are six tau isoforms which are formed by alternate splicing. Hyperphosphorylation of Tau protein results in aggregation and formation of neurofibrillary tangles, which is significantly associated with cognitive impairment in Alzheimer’s disease. Tau Thr181 phosphorylation has been used as a biomarker for the differentiation of Alzheimer’s disease and other primary dementia such as dementia with Lewy bodies.

Applications: **Indirect ELISA:** The antibody reacts with phosphorylated peptide [Ala-Pro-Lys-Thr(p)-Pro-Pro-Ser-Ser-Gly-Glu-Cys] coated ELISA plate in indirect ELISA.

Western Blot: The antibody detected a 46kD band in mouse brain tissue lysate.

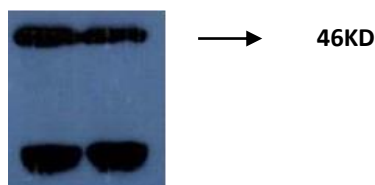

The antibody (Tau-pho181 9B4) was used at 1:2000 (0.5 µg/mL) dilution for Western Blot.

References: If research is published using this product, please inform Anogen in order to cite the reference on this datasheet. Anogen will provide one unit of product in the same category as gratitude.

This product is for LABORATORY RESEARCH USE and further manufacture ONLY, and cannot be administrated to human and animals for use in diagnostic and therapeutic procedures.

Manufactured by ANOGEN - A Division of YES Biotech Laboratories Ltd.